




OWANDY-CAM^{HD}

ENJOY EASIER COMMUNICATION WITH YOUR PATIENTS
DIGITAL USB INTRAORAL CAMERA

Owandy-Cam

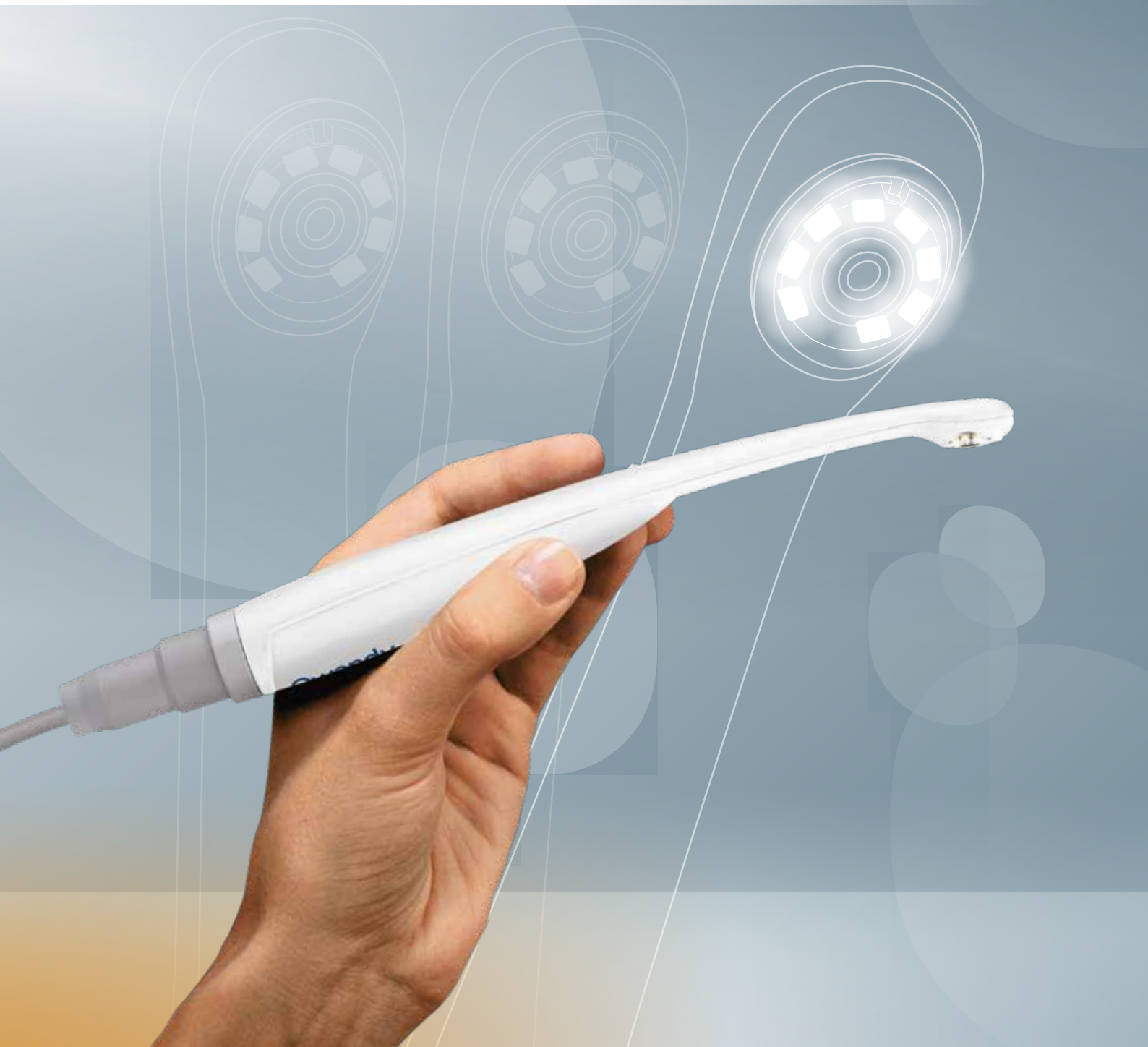

Facilitate communication with patients and treatment plan acceptance

Explain clinical procedures to your patients more easily!
Your patients will feel more involved and have greater trust in the process.



4 shot modes


Its "focus free" feature gives a greater field depth and high quality images. Requiring no adjustments at all, this unit delivers perfectly clear images, from a full smile to just one tooth - and will save you time, too!

1



SMILE

2

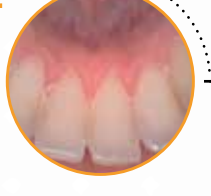

INTRAORAL

3



INTRAORAL

4



MACRO



100% integrated
The Owandy-Cam HD is fully embedded in all leading software packages.





Ergonomic design and wide angle de 105° lens

Its curved shape and ultra-slim head make for greater patient comfort and prevent injuries to the mouth.

Its 105° viewing angle sets it apart from other cameras on the market, by enabling more effective exploration of distal areas.



User-friendly

Set your images by simply sliding over the touch button. There'll be no blurry images, allowing you to work with complete confidence.

The USB2 connection makes it easy to move the camera around and share it with your colleagues and partners, anywhere in your surgery!

Technical features

Sensor	High sensitivity 1/4 " Sony CCD
Resolution	50 images / sec
Lighting	8 LEDs
Adjustment	Focus free, large field depth
View angle	105°
USB cable length	2.5m
Handpiece size (L x W x H)	205 x 28 x 24mm
Thickness of intraoral part	11mm
Handset weight	55g
Image setter	Touch button



"Push&Pull" system

The Owandy-Cam HD fits easily onto your dental chair, just like all your handpieces.



DIGITAL WORKFLOW OWANDY RADIOLOGY

A COMPREHENSIVE RANGE TO MEET ALL YOUR REQUIREMENTS

